Product information

Information about other products is available at: www.demeditec.com

Nephelines Urine ELISA

DEE8600

2 x 96

Demeditec Diagnostics GmbH
Lise-Meitner-Strasse 2
24145 Kiel – Germany
www.demeditec.com

Content / Inhaltsverzeichnis

1. Introduction	3
2. Procedural cautions, guidelines, warnings and limitations	3
3. Storage and stability	4
4. Materials	4
5. Sample collection and storage	6
6. Test procedure	6
7. Calculation of results	8
8. Assay characteristics	9
9. References/Literature	10

SYMBOLS USED WITH DEMEDITEC ASSAYS 10

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1. Introduction

1.1 Intended use and principle of the test

Enzyme Immunoassay for the quantitative determination of Metanephrine and Normetanephrine in urine.

During the sample preparation Metanephrine (Metadrenaline) and Normetanephrine (Normetadrenaline) are quantitatively acylated.

The subsequent competitive ELISA kit uses the microtiter plate format. The antigen is bound to the solid phase of the microtiter plate. The acylated standards, controls and samples and the solid phase bound analytes compete for a fixed number of antibody binding sites. After the system is in equilibrium, free antigen and free antigen-antibody complexes are removed by washing. The antibody bound to the solid phase is detected by an anti-rabbit IgG-peroxidase conjugate using TMB as a substrate. The reaction is monitored at 450 nm.

Quantification of unknown samples is achieved by comparing their absorbance with a reference curve prepared with known standard concentrations.

⚠ *The anti-Metanephrine antibodies used in this test kit only recognise the biologically relevant L-forms of Metanephrine. Commercially available synthetic Metanephrine is always a mixture of the D- and L-forms. The ratio between both forms differs widely from lot to lot. This has important implications if synthetic Metanephrine is used to enrich native samples. As only about 50% of the synthetic Metanephrine – the L-portion – will be detected by use of this kit, spiked samples will be underestimated. Therefore, native samples containing solely the L-form should be used.*

1.2 Clinical application

Metanephrine and Normetanephrine are the metabolites of the catecholamines Epinephrine and Norepinephrine, respectively. They are metabolized to Vanillylmandelic acid or excreted with the urine. Patients with pheochromocytoma or other tumors derived from neuroendocrine cells show elevated urinary levels of total Metanephrines.

As catecholamine secretion from neuroendocrine cells might show high variations, urine samples collected over a period of 24 hours are used to average these fluctuations.

Therapeutic consequences should never be based on laboratory results alone even if all test results are in agreement with the items as under point "Procedural cautions, guidelines and warnings". Any laboratory result is only a part of the total clinical picture of the patient.

Only in cases where the laboratory results are in an acceptable agreement with the overall clinical picture of the patient it can be used for therapeutic consequences.

The test result itself should never be the sole determinant for deriving any therapeutic consequences.

2. Procedural cautions, guidelines, warnings and limitations

2.1 Procedural cautions, guidelines and warnings

- (1) This kit is intended for professional use only. Users should have a thorough understanding of this protocol for the successful use of this kit. Only the test instruction provided with the kit is valid and has to be used to run the assay. Reliable performance will only be attained by strict and careful adherence to the instructions provided.
- (2) This assay was validated for a certain type of sample as indicated in *Intended Use* (please refer to Chapter 1). Any off-label use of this kit is in the responsibility of the user and the manufacturer cannot be held liable.
- (3) The principles of Good Laboratory Practice (GLP) have to be followed.
- (4) In order to reduce exposure to potentially harmful substances, wear lab coats, disposable protective gloves and protective glasses where necessary.
- (5) All kit reagents and specimens should be brought to room temperature and mixed gently but thoroughly before use. Avoid repeated freezing and thawing of reagents and specimens.
- (6) For dilution or reconstitution purposes, use deionized, distilled, or ultra-pure water.
- (7) The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch with desiccant and used in the frame provided. Microtiter strips which are removed from the frame for usage should be marked accordingly to avoid any mix-up.
- (8) Duplicate determination of sample is highly recommended to be able to identify potential pipetting errors.
- (9) Once the test has been started, all steps should be completed without interruption. Make sure that the required reagents, materials and devices are prepared ready at the appropriate time.
- (10) Incubation times do influence the results. All wells should be handled in the same order and time intervals.

- (11) To avoid cross-contamination of reagents, use new disposable pipette tips for dispensing each reagent, sample, standard and control.
- (12) A standard curve must be established for each run.
- (13) The controls should be included in each run and fall within established confidence limits. The confidence limits are listed in the QC-Report.
- (14) Do not mix kit components with different lot numbers within a test and do not use reagents beyond expiry date as shown on the kit labels.
- (15) Avoid contact with Stop Solution containing 0.25 M H₂SO₄. It may cause skin irritation and burns. In case of contact with eyes or skin, rinse off immediately with water.
- (16) TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them.
- (17) For information on hazardous substances included in the kit please refer to Safety Data Sheet (SDS). The Safety Data Sheet for this product is made available directly on the website of the manufacturer or upon request.
- (18) The expected reference values reported in this test instruction are only indicative. It is recommended that each laboratory establishes its own reference intervals.
- (19) The results obtained with this test kit should not be taken as the sole reason for any therapeutic consequence (e.g. medication before a scheduled surgery) but have to be correlated to other diagnostic tests and clinical observations.
- (20) Kit reagents must be regarded as hazardous waste and disposed according to national regulations.
- (21) In case of any severe damage to the test kit or components, the manufacturer has to be informed in writing, at the latest, one week after receiving the kit. Severely damaged single components must not be used for a test run. They must be stored properly until the manufacturer decides what to do with them. If it is decided that they are no longer suitable for measurements, they must be disposed of in accordance with national regulations.

2.2 Limitations

Any inappropriate handling of samples or modification of this test might influence the results.

2.2.1 Interfering substances

24-hour urine

Please note the sample preparation! If the percentage of the final concentration of acid is too high, this will lead to incorrect results for the urine samples.

2.2.2 Drug interferences

There are no known substances (drugs) which ingestion interferes with the measurement of (Nor-)metanephrine level in the sample.

2.2.3 High-Dose-Hook effect

No hook effect was observed in this test.

3. Storage and stability

Store the unopened reagents at 2–8 °C until expiration date. Do not use components beyond the expiry date indicated on the kit labels. Once opened the reagents are stable for 2 months when stored at 2–8 °C. Once the resealable pouch has been opened, care should be taken to close it tightly with desiccant again.

4. Materials

4.1 Content of the kit

BA D-0023	REAC-TUBES	Reaction Tubes – Ready to use
Content:	Reaction Tubes in a resealable pouch	
Volume:	2 x 50 tubes	
BA E-0030	WASH-CONC 50x	Wash Buffer Concentrate – Concentrated 50x
Content:	Buffer with a non-ionic detergent and physiological pH	
Volume:	1 x 20 ml/vial, light purple cap	

- BA E-0045** **CONJUGATE** **Enzyme Conjugate** – Ready to use
 Content: Goat anti-rabbit immunoglobulins conjugated with peroxidase
 Volume: 2 x 12 ml/vial, red cap
- BA E-0055** **SUBSTRATE** **Substrate** – Ready to use
 Content: Chromogenic substrate containing tetramethylbenzidine, substrate buffer and hydrogen peroxide
 Volume: 2 x 12 ml/vial, black cap
- BA E-0080** **STOP-SOLN** **Stop Solution** – Ready to use
 Content: 0.25 M sulfuric acid
 Volume: 2 x 12 ml/vial, light grey cap
- Hazards identification: 

H290 May be corrosive to metals.

- BA E-0131** **ADR MN** **Metanephrine Microtiter Strips** – Ready to use
 Content: 1 x 96 well (12x8) antigen precoated microwell plate in a resealable blue pouch with desiccant
- BA E-0231** **NAD NMN** **Normetanephrine Microtiter Strips** – Ready to use
 Content: 1 x 96 well (12x8) antigen precoated microwell plate in a resealable yellow pouch with desiccant
- BA E-8410** **MN-AS** **Metanephrine Antiserum** – Ready to use
 Content: Rabbit anti-Metanephrine antibody, blue coloured
 Volume: 1 x 12 ml/vial, blue cap
- BA E-8510** **NMN-AS** **Normetanephrine Antiserum** – Ready to use
 Content: Rabbit anti-Normetanephrine antibody, yellow coloured
 Volume: 1 x 12 ml/vial, yellow cap

Standards and Controls – Ready to use

Cat. no.	Component	Colour/Cap	Concentration ng/ml (= µg/l)		Concentration nmol/l		Volume/ Vial
			MN	NMN	MN	NMN	
BA R-8601	STANDARD A	white	0	0	0	0	4 ml
BA R-8602	STANDARD B	light yellow	20	30	101	164	4 ml
BA R-8603	STANDARD C	orange	60	90	304	491	4 ml
BA R-8604	STANDARD D	dark blue	200	300	1014	1638	4 ml
BA R-8605	STANDARD E	light grey	600	900	3042	4914	4 ml
BA R-8606	STANDARD F	black	2000	3000	10140	16380	4 ml
BA R-8651	CONTROL 1	light green	Refer to QC-Report for expected value and acceptable range!				4 ml
BA R-8652	CONTROL 2	dark red	Refer to QC-Report for expected value and acceptable range!				4 ml

Conversion: Metanephrine (ng/ml) x 5.07 = Metanephrine (nmol/l)
 Normetanephrine (ng/ml) x 5.46 = Normetanephrine (nmol/l)

Content: Acidic buffer with non-mercury preservatives, spiked with defined quantity of Metanephrine and Normetanephrine

- BA R-0012** **ACYL-CONC** **Acylation Concentrate** – Concentrated
 Content: Concentrated acylation reagent
 Volume: 1 x 0.5 ml/vial, white cap

Hazards identification: 

H314 Causes severe skin burns and eye damage.

BA R-0075	ACYL-DILUENT	Acylation Diluent – Ready to use
Content:	Dimethylsulfoxide	
Volume:	1 x 3.5 ml/vial, white cap	
BA R-8611	ACYL-BUFF	Acylation Buffer – Ready to use
Content:	TRIS buffer	
Volume:	1 x 30 ml/vial, white cap	
BA R-8619	HCL	Hydrochloric Acid – Ready to use
Content:	0.25 M hydrochloric acid, yellow coloured	
Volume:	1 x 30 ml/vial, dark green cap	

4.2 Additional materials and equipment required but not provided in the kit

- Calibrated precision pipettes to dispense volumes between 10–600 µl; 1.2–3 ml
- Microtiter plate washing device (manual, semi-automated or automated)
- ELISA reader capable of reading absorbance at 450 nm and if possible 620–650 nm
- Absorbent material (paper towel)
- Water (deionized, distilled, or ultra-pure)
- Vortex mixer
- Temperature controlled water bath (90 °C) or similar heating device

⚠ The assay can be performed with or without shaking. If a microtiter plate shaker is used, it should have the following characteristics: shaking amplitude 3 mm; approx. 600 rpm

5. Sample collection and storage

Spontaneous or 24-hour urine, collected in a bottle containing 10–15 ml of 6 M HCl, should be used. Determine the total volume of urine excreted during a period of 24 h for calculation of the results. Storage: up to 5 days at 2–8 °C, for longer periods (up to 6 months) at -20 °C. Repeated freezing and thawing should be avoided. Avoid exposure to direct sunlight.

6. Test procedure

Allow all reagents to reach room temperature and mix thoroughly by gentle inversion before use. Number the Reaction Tubes accordingly. Duplicate determinations are recommended. It is recommended to number the strips of the microwell plate before usage to avoid any mix-up.

⚠ The sample preparation (hydrolysis and acylation) is identical for both the Metanephrine and Normetanephrine assay and has to be done only once.

The binding of the antisera and of the enzyme conjugate and the activity of the enzyme are temperature dependent. The higher the temperature, the higher the absorption values will be. Varying incubation times will have similar influences on the absorbance. The optimal temperature during the Enzyme Immunoassay is between 20–25 °C.

6.1 Preparation of reagents

Wash Buffer

Dilute the 20 ml Wash Buffer Concentrate with water (deionized, distilled, or ultra-pure) to a final volume of 1000 ml. Storage: 2 months 2–8 °C

Acylation Solution

⚠ Before preparing the Acylation Solution make sure that the Acylation Diluent (BA R-0075) has reached room temperature (≥ 20 °C) and forms a homogenous, crystal-free solution.

Dilute the Acylation Concentrate (BA R-0012) 1 + 60 with Acylation-Diluent in a glass or polypropylene-vial.

Acylation Concentrate	10 µl	15 µl	25 µl	50 µl
Acylation-Diluent	600 µl	900 µl	1.5 ml	3 ml

⚠ The Acylation Solution has to be prepared freshly prior to the assay (not longer than 60 minutes in advance). Discard after use!

Metanephrine Microtiter Strips and Normetanephrine Microtiter Strips

In rare cases residues of the blocking and stabilizing reagent can be seen in the wells as small, white dots or lines. These residues do not influence the quality of the product.

6.2 Sample preparation and acylation**Hydrolysis**

- | | |
|----|---|
| 1. | Pipette 25 µl of standards, controls, and urine samples into the respective Reaction Tubes . |
| 2. | Add 250 µl Hydrochloric Acid to all tubes. |
| 3. | Mix thoroughly (vortex) and hydrolyze for 30 min at 90 °C . |
| 4. | Cool down the tubes to room temperature. |

Acylation

- | | |
|---|---|
| 1. | Pipette 250 µl of Acylation Buffer into all tubes. |
| 2. | Add 25 µl of Acylation Solution (<i>refer to 6.1</i>) to all tubes. |
| 3. | Mix thoroughly (vortex) and acylate for 15 min at RT (20–25 °C) . |
| 4. | Add 2.5 ml water (deionized, distilled, or ultra-pure) to all tubes. |
|  | Take 25 µl of the acylated standards, controls and urine samples for the Metanephrine ELISA and Normetanephrine ELISA . |

6.3 Metanephrine ELISA

The usage of a shaker is not mandatory. The alternative protocol without shaker is highlighted in italic and shaded in grey.

- | | |
|-----|--|
| 1. | Pipette 25 µl of the acylated standards, controls and samples into the appropriate wells of the Metanephrine Microtiter Strips . |
| 2. | Pipette 100 µl of the Metanephrine Antiserum into all wells. |
| 3. | Incubate 30 min at RT (20–25 °C) on a shaker (approx. 600 rpm).
<i>Without usage of a shaker: shake the Metanephrine Microtiter Strips shortly by hand and incubate for 1 h at RT (20–25 °C).</i> |
| 4. | Discard or aspirate the content of the wells. Wash the plate 3 x by adding 300 µl of Wash Buffer , discarding the content and blotting dry each time by tapping the inverted plate on absorbent material. |
| 5. | Pipette 100 µl of the Enzyme Conjugate into all wells. |
| 6. | Incubate for 15 min at RT (20–25 °C) on a shaker (approx. 600 rpm).
<i>Without usage of a shaker: incubate for 15 min at RT (20–25 °C).</i> |
| 7. | Discard or aspirate the content of the wells. Wash the plate 3 x by adding 300 µl of Wash Buffer , discarding the content and blotting dry each time by tapping the inverted plate on absorbent material. |
| 8. | Pipette 100 µl of the Substrate into all wells. |
| 9. | Incubate for 15 ± 2 min at RT (20–25 °C) on a shaker (approx. 600 rpm).
<i>Without usage of a shaker: incubate for 15 min ± 2 at RT (20–25 °C).</i>
 Avoid exposure to direct sunlight! |
| 10. | Add 100 µl of the Stop Solution to each well and shake the microtiter plate to ensure a homogeneous distribution of the solution. |
| 11. | Read the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to 450 nm (if available a reference wavelength between 620 nm and 650 nm is recommended). |

6.4 Normetanephrine ELISA

The usage of a shaker is not mandatory. The alternative protocol without shaker is highlighted in italic and shaded in grey.

1.	Pipette 25 µl of the acylated standards, controls and samples into the appropriate wells of the Normetanephrine Microtiter Strips .
2.	Pipette 100 µl of the Normetanephrine Antiserum into all wells.
3.	Incubate 30 min at RT (20–25 °C) on a shaker (approx. 600 rpm). <i>Without usage of a shaker: shake the Normetanephrine Microtiter Strips shortly by hand and incubate for 1 h at RT (20–25 °C).</i>
4.	Discard or aspirate the content of the wells. Wash the plate 3 x by adding 300 µl of Wash Buffer , discarding the content and blotting dry each time by tapping the inverted plate on absorbent material.
5.	Pipette 100 µl of the Enzyme Conjugate into all wells.
6.	Incubate for 15 min at RT (20–25 °C) on a shaker (approx. 600 rpm). <i>Without usage of a shaker: incubate for 15 min at RT (20–25 °C).</i>
7.	Discard or aspirate the content of the wells. Wash the plate 3 x by adding 300 µl of Wash Buffer , discarding the content and blotting dry each time by tapping the inverted plate on absorbent material.
8.	Pipette 100 µl of the Substrate into all wells.
9.	Incubate for 15 ± 2 min at RT (20–25 °C) on a shaker (approx. 600 rpm). <i>Without usage of a shaker: incubate for 15 min ± 2 at RT (20–25 °C).</i>
▲	Avoid exposure to direct sun light!
10.	Add 100 µl of the Stop Solution to each well and shake the microtiter plate to ensure a homogeneous distribution of the solution.
11.	Read the absorbance of the solution in the wells within 10 minutes, using a microplate reader set to 450 nm (if available a reference wavelength between 620 nm and 650 nm is

7. Calculation of results

Measuring range	Metanephrine	Normetanephrine
	10.5–2000 ng/ml	16.2–3000 ng/ml

The standard curve is obtained by plotting the absorbance readings (calculate the mean absorbance) of the standards (linear, y-axis) against the corresponding standard concentrations (logarithmic, x-axis). Use a non-linear regression for curve fitting (e.g. 4-parameter, marquardt).

⚠ *This assay is a competitive assay. This means: the OD-values are decreasing with increasing concentrations of the analyte. OD-values found below the standard curve correspond to high concentrations of the analyte in the sample and have to be reported as being positive.*

The concentrations of the samples and controls can be read directly from the standard curve.

The amount of analyte excreted per day (µg/day) is calculated according to:

concentration of the sample (in µg/l) x volume of urine excreted per day (in l/day)

Example

The concentration of the sample read from the curve is 125 µg/l. The amount of urine collected during 24 hours is 1.3 l. Then the amount of analyte excreted during one day would be:

$$125 \text{ µg/l} \times 1.3 \text{ l/day} = 162.5 \text{ µg/day}$$

Samples found with concentrations higher than the highest standard (Standard F) should be diluted accordingly with Standard A and have to be re-assayed.

Conversion

Metanephrine (ng/ml) x 5.07 = Metanephrine (nmol/l)

Normetanephrine (ng/ml) x 5.46 = Normetanephrine (nmol/l)

Expected reference values

It is strongly recommended that each laboratory should determine its own reference values.

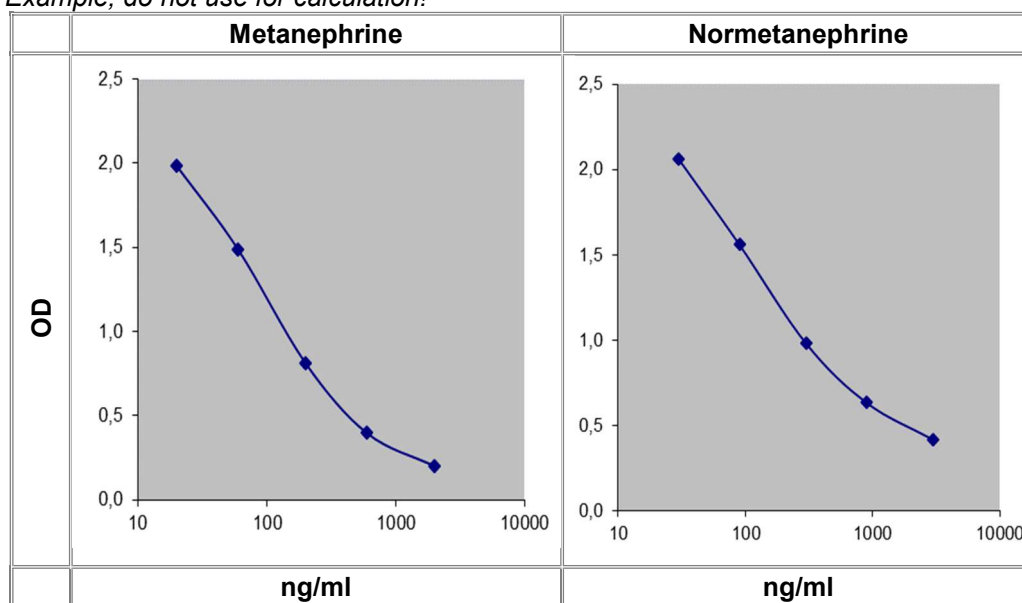
	Metanephrine	Normetanephrine
24-hour urine	< 350 µg/day	< 600 µg/day

7.1 Quality control

It is recommended to use control samples according to national regulations. Use controls at both normal and pathological levels. The kit or other commercial controls should fall within established confidence limits. The confidence limits of the kit controls are indicated on the QC-Report.

7.2 Typical standard curves

⚠ *Example, do not use for calculation!*

**8. Assay characteristics**

Analytical Sensitivity		Metanephrine	Normetanephrine
	LOD	8.6 ng/ml	14.7 ng/ml
	LOB	6.0 ng/ml	10.4 ng/ml
	LOQ	10.5 ng/ml	16.2 ng/ml

Analytical Specificity (Cross Reactivity)	Substance	Cross Reactivity (%)	
		Metanephrine	Normetanephrine
	Derivatized Metanephrine	100	0.11
	Derivatized Normetanephrine	0.15	100
	Derivatized 3-methoxytyramine	< 0.01	0.19
	Adrenaline	3.3	< 0.01
	Noradrenaline	< 0.01	0.64
	Dopamine	< 0.01	< 0.01
	Vanillic mandelic acid, L-Dopa, Homovanillic acid, L-Tyrosin, Tyramin	< 0.01	< 0.01

Precision							
Intra-Assay				Inter-Assay			
	Sample	Range (ng/ml)	CV (%)		Sample	Range (ng/ml)	CV (%)
Metanephrine	1	34.4 ± 3.1	9	Metanephrine	1	32.8 ± 5.4	16
	2	59.8 ± 5.1	9		2	57.2 ± 9.0	16
	3	141 ± 13.8	10		3	144 ± 25.0	17
	4	575 ± 71.4	12		4	394 ± 64.1	16

Normetanephrine	1	53.9 ± 7.1	13	Normetanephrine	1	52.5 ± 10.7	20
	2	116 ± 12.3	11		2	103 ± 15.2	15
	3	322 ± 28.1	9		3	276 ± 38.7	14
	4	1121 ± 128	11		4	738 ± 83.0	11

Linearity		Serial dilution up to	Mean (%)	Range (%)
	Metanephrine	1:64	102	94–115
	Normetanephrine	1:64	101	90–113

Recovery		Mean (%)	Range (%)	Range (ng/ml)
	Metanephrine	97	85–113	20.2–1484
	Normetanephrine	100	93–111	26.5–3124

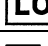
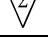
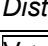
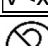
Method Comparison versus HPLC*	Metanephrine	HPLC = 0.9 ELISA – 0.8	r = 0.99; n = 40
	Normetanephrine	HPLC = 0.9 ELISA + 0.6	r = 0.99; n = 40
*The concentrations were assessed using both the ELISA and the HPLC method (external QC samples from UK NEQAS). The correlation between ELISA and HPLC is excellent. Please take in mind, that the UK control values are the mean of about 40 different HPLC users, and contain always one pathological sample per sending.			

9. References/Literature

- (1) Parrott et al. Urinary corticosterone and normetanephrine levels after voluntary wheel and forced treadmill running in the db/db mouse. *Journal of Diabetes Mellitus*, 1(4):71-78 (2011)
- (2) Petramala et al. Multiple Catecholamine-Secreting Paragangliomas: Diagnosis after Hemorrhagic Stroke in a Young Woman. *Endocrine Practice*, 14(3):340-346 (2008)
- (3) Sato et al. Central control of bone remodeling by neuromedin U. *Nature Medicine*, 13:1234-1240 (2007)

For updated literature or any other information please contact your local supplier.

SYMBOLS USED WITH DEMEDITEC ASSAYS

Symbol	English	Deutsch	Française	Espanol	Italiano
	European Conformity	CE-Konformitätskennzeichnung	Conforme aux normes européennes	Conformidad europea	Conformità europea
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las Instrucciones	Consultare le istruzioni per l'uso
	In vitro diagnostic device	In-vitro-Diagnostikum	utilisation Diagnostic in vitro	Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Référence	Número de catálogo	No. di catalogo
	Lot. No. / Batch code	Chargen-Nr.	No. de lot	Número de lote	Lotto no
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Note warnings and precautions	Warnhinweise und Vorsichtsmaßnahmen beachten	Avertissements et mesures de précaution font attention	Tiene en cuenta y advertencias precauciones	Annoti avvisi e le precauzioni
	Storage Temperature	Lagerungstemperatur	Température de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeitsdatum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
<i>Distributed by</i>	Distributed by	Vertrieb durch	Distribution par	Distribución por	Distribuzione da parte di
	Version	Version	Version	Versión	Versione
	Single-use	Einmalverwendung	À usage unique	Uso único	Uso una volta

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