# Product information

Information about other products is available at: [www.demeditec.com](http://www.demeditec.com)

User's Manual

# AMH ELISA

**DE6141**

**96 wells**

Demeditec Diagnostics GmbH  
Lise-Meitner-Strasse 2  
24145 Kiel – Germany  
[www.demeditec.com](http://www.demeditec.com)

**Please use only the valid version of the Instructions for Use provided with the kit.  
 Verwenden Sie nur die jeweils gültige, im Testkit enthaltene, Gebrauchsanweisung.  
 Si prega di usare la versione valida delle istruzioni per l'uso a disposizione con il kit.  
 Por favor, use sólo la versión válida de las instrucciones de uso que se suministran con el kit.  
 Utilisez seulement la version valide des Instructions d'utilisation fournies avec le kit.**

**Table of Contents / Inhaltsverzeichnis / Tabella dei Contenuti / Tabla de Contenidos / Table des matières**

1	INTENDED USE .....	4
2	PRINCIPLE OF THE TEST .....	4
3	WARNINGS AND PRECAUTIONS .....	5
4	REAGENTS.....	6
5	SPECIMEN COLLECTION AND PREPARATION .....	7
6	ASSAY PROCEDURE .....	7
7	EXPECTED NORMAL VALUES .....	9
8	QUALITY CONTROL .....	9
9	PERFORMANCE CHARACTERISTICS .....	10
10	LIMITATIONS OF USE .....	12
11	LEGAL ASPECTS .....	12
1	ZWECKBESTIMMUNG .....	13
2	TESTPRINZIP .....	13
3	VORSICHTSMAßNAHMEN .....	13
4	BESTANDTEILE DES KITS .....	14
5	PROBENVORBEREITUNG .....	15
6	TESTDURCHFÜHRUNG .....	15
7	ERWARTETE WERTE.....	17
8	QUALITÄTSKONTROLLE.....	17
9	ASSAY-CHARAKTERISTIKA.....	18
10	GRENZEN DES TESTS.....	18
11	RECHTLICHE GRUNDLAGEN .....	19
1	DESTINAZIONE D'USO .....	20
2	PRINCIPIO DEL TEST .....	20
3	PRECAUZIONI .....	20
4	COMPONENTI DEL KIT .....	21
5	CAMPIONI.....	21
6	ATTUAZIONE DEL TEST .....	22
7	VALORI NORMALI.....	24
8	CONTROLLO QUALITÀ.....	24
9	CARATTERISTICHE DEL TEST.....	24
10	LIMITAZIONE DEL TEST .....	25
11	ASPETTI LEGALI .....	25

---

1	FINALIDAD PREVISTA.....	26
2	FUNDAMENTO DEL ENSAYO .....	26
3	PRECAUCIONES.....	26
4	COMPONENTES DEL KIT.....	27
5	MUESTRAS.....	27
6	PROCEDIMIENTO DE ENSAYO .....	28
7	VALORES ESPERADOS .....	30
8	CONTROL DE CALIDAD .....	30
9	CARACTERÍSTICAS DEL ENSAYO .....	31
10	LIMITACIONES DE USO .....	31
11	ASPECTOS LEGALES .....	32
1	UTILISATION PRÉVUE .....	33
2	PRINCIPE DU TEST .....	33
3	AVERTISSEMENTS ET PRECAUTIONS.....	33
4	RÉACTIFS.....	34
5	PRÉLÈVEMENT ET PRÉPARATION DES SPÉCIMENS .....	34
6	RÉALISATION DU TEST .....	36
7	VALEURS NORMALES ATTENDUES .....	37
8	CONTRÔLE DE QUALITÉ .....	37
9	CARACTERISTIQUES DES PERFORMANCES.....	37
10	LIMITES D'UTILISATION.....	38
11	ASPECTS JURIDIQUES.....	39
12	REFERENCES / LITERATURE .....	40
	SYMBOLS USED WITH DEMEDITEC ASSAYS .....	40

## 1 INTENDED USE

The **DEMEDIATEC AMH ELISA** is an enzyme immunoassay for the quantitative **in vitro diagnostic** measurement of Anti-Müllerian Hormone (AMH) in serum or plasma (EDTA or lithium heparin plasma)

### 1.1 Summary and Explanation

Anti-Müllerian Hormone (AMH), a dimeric 140 KDa glycoprotein, is a member of the transforming growth factor- $\beta$  (TGF- $\beta$ ) family of cytokines which plays an essential role in the normal differentiation of reproductive structures. (1) AMH is secreted by Sertoli cells of the testes during embryogenesis of the fetal male, preventing the development of the Müllerian ducts to the uterus and other Müllerian structures. In females, AMH is secreted by the granulosa cells of ovarian follicles. (2,3) AMH has been identified as a reliable marker of ovarian reserve that can help predict early ovarian follicle loss and menopause onset. AMH levels also reflect the effects of damaging gynecologic surgeries or gonadotoxic treatments such as chemotherapy on ovarian reserve. Furthermore, AMH contributes in the diagnosis of certain diseases such as granulosa cell tumors or polycystic ovary syndrome (PCOS). (4,5,6,7,8,9,10)

## 2 PRINCIPLE OF THE TEST

The DEMEDIATEC AMH ELISA is a solid phase enzyme-linked immunosorbent assay (ELISA) based on the **sandwich principle**. The microtiter wells are coated with a monoclonal (mouse) antibody directed towards a unique antigenic site of the AMH molecule. During the first incubation, the AMH in the added sample bind to the immobilized antibody. The simultaneously added enzyme conjugate, which contains an AMH antibody conjugated to horseradish peroxidase, binds to the AMH forming a sandwich complex. After a washing step to remove all unbound substances, the solid phase is incubated with the substrate solution. The colorimetric reaction is stopped by addition of stop solution, and optical density (OD) of the resulting yellow product is measured. The intensity of colour is proportional to the concentration of the analyte in the sample. A standard curve is constructed by plotting OD values against concentrations of standards, and concentrations of unknown samples are determined using this standard curve.

### 3 WARNINGS AND PRECAUTIONS

1. This kit is for in vitro diagnostic use only. For professional use only.
2. All reagents of this test kit which contain human serum or plasma have been tested and confirmed negative for HIV I/II, HBsAg and HCV by FDA approved procedures. All reagents, however, should be treated as potential biohazards in use and for disposal.
3. Before starting the assay, read the instructions completely and carefully. Use the valid version of instructions for use provided with the kit. Be sure that everything is understood.
4. The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch and used in the frame provided.
5. Pipetting of samples and reagents must be done as quickly as possible and in the same sequence for each step.
6. Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution coloured. Do not pour reagents back into vials as reagent contamination may occur.
7. Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells.
8. Do not let wells dry during assay; add reagents immediately after completing the rinsing steps.
9. Allow the reagents to reach room temperature (20 °C to 26 °C) before starting the test. Temperature will affect the optical density readings of the assay. However, values for the patient samples will not be affected.
10. Never pipet by mouth and avoid contact of reagents and specimens with skin and mucous membranes.
11. Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.
12. Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents or specimens may give false results.
13. Handling should be done in accordance with the procedures defined by an appropriate national biohazard safety guideline or regulation.
14. Do not use reagents beyond expiry date as shown on the kit labels.
15. All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microtiter plate readers.
16. Do not mix or use components from kits with different lot numbers. It is advised not to exchange wells of different plates even of the same lot. The kits may have been shipped or stored under different conditions and the binding characteristics of the plates may result slightly different.
17. Avoid contact with Stop Solution containing 0.5 M H<sub>2</sub>SO<sub>4</sub>. It may cause skin irritation and burns.
18. Some reagents contain Proclin 300, BND and/or MIT as preservatives. In case of contact with eyes or skin, flush immediately with water.
19. TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them. If inhaled, take the person to open air.
20. Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guideline or regulation.
21. For information on hazardous substances included in the kit please refer to Safety Data Sheets. Safety Data Sheets for this product are available upon request directly from DEMEDITEC Diagnostics GmbH.

## 4 REAGENTS

### 4.1 Reagents provided

1. **SORB MT Microtiterwells**, 12 x 8 (break apart) strips, 96 wells;  
Wells coated with anti-AMH antibody (monoclonal).
2. **CAL 0 - 5 Standard (Standard 0 - 5)**, 6 vials, 1 mL each, ready to use;  
Concentrations: 0.0 – 0.4 – 1.0 – 4.0 – 10 – 20 ng/mL  
Conversion: 1 ng/mL = 7.14 pmol/L  
The standards are calibrated against the following reference material: 1<sup>st</sup> WHO International Reference Reagent, Mullerian Inhibiting Substance/Anti-Mullerian Hormone, NIBSC code: 16/190  
Contain non-mercury preservative.
3. **CONTROL low & high Control Low & High**, 2 vials, 1 mL each, ready to use; For control values and ranges please refer to Certificate of Analysis. Contain non-mercury preservative.
4. **ENZ CONJ Enzyme Conjugate**, 1 vial, 14 mL, ready to use; Anti-AMH antibody conjugated with horseradish peroxidase; Contains non-mercury preservative.
5. **SUB TMB Substrate Solution**, 1 vial, 14 mL, ready to use; Tetramethylbenzidine (TMB).
6. **STOP SOLN Stop Solution**, 1 vial, 14 mL, ready to use; Contains 0.5 M H<sub>2</sub>SO<sub>4</sub>. Avoid contact with the stop solution. It may cause skin irritations and burns.
7. **WASH SOLN 40x Wash Solution**, 1 vial, 30 mL (40X concentrated); See "Reagent Preparation".

**Note:** Additional Standard 0 for sample dilution is available upon request.

### 4.2 Materials required but not provided

- A calibrated microtiter plate reader (450 nm, with reference wavelength at 620 nm to 630 nm)
- Calibrated variable precision micropipettes
- Absorbent paper
- Distilled water
- Timer
- Graph paper or software for data reduction

### 4.3 Storage Conditions

When stored at 2 °C to 8 °C unopened reagents will retain reactivity until expiration date. Do not use reagents beyond this date. Opened reagents must be stored at 2 °C to 8 °C. Microtiter wells must be stored at 2 °C to 8 °C. Once the foil bag has been opened, care should be taken to close it tightly again. Opened kits retain activity for 8 weeks if stored as described above.

### 4.4 Reagent Preparation

Bring all reagents and required number of strips to room temperature (20 °C to 26 °C) prior to use.

### Wash Solution

Add distilled water to the 40X concentrated Wash Solution. Dilute 30 mL of concentrated Wash Solution with 1170 mL distilled water to a final volume of 1200 mL. The diluted Wash Solution is stable for 2 weeks at room temperature.

### 4.5 Disposal of the Kit

The disposal of the kit and all used materials/reagents must be performed according to the national regulations. Special information for this product is given in the Safety Data Sheet, section 13.

### 4.6 Damaged Test Kits

In case of any damage to the test kit or components, DEMEDITEC must be informed in writing, at the latest one week after receiving the kit. Damaged single components should not be used for a test run. They have to be stored until a final solution has been found. After this, they should be disposed of according to the official regulations.

## 5 SPECIMEN COLLECTION AND PREPARATION

Serum or plasma (EDTA or lithium heparin) can be used in this assay.

**Note:** Samples containing sodium azide should not be used in the assay.

In general, it should be avoided to use haemolytic, icteric, or lipaemic specimens. For further information refer to chapter "Interfering Substances".

### 5.1 Specimen Collection

#### **Serum:**

Collect blood by venipuncture (e.g. Sarstedt Monovette for serum), allow to clot, and separate serum by centrifugation at room temperature. Do not centrifuge before complete clotting has occurred. Patients receiving anticoagulant therapy may require increased clotting time.

#### **Plasma:**

Whole blood should be collected into centrifuge tubes containing anticoagulant (e.g. Sarstedt Monovette with the appropriate plasma preparation) and centrifuged immediately after collection.

### 5.2 Specimen Storage and Preparation

Specimens should be capped and may be stored for up to 7 days at 2 °C to 8 °C prior to assaying.

Specimens stored for a longer time (up to 3 months) should be frozen only once at -20 °C prior to assay.

Thawed samples should be inverted several times prior to testing.

### 5.3 Specimen Dilution

If in an initial assay, a specimen is found to contain more analyte than the highest standard, the specimen can be diluted with Standard 0 and re-assayed as described in "Assay Procedure".

For the calculation of the concentrations this dilution factor has to be taken into account.

#### Example:

a) dilution 1:10:           10 µL sample + 90 µL Standard 0 (mix thoroughly)

b) dilution 1:100:        10 µL dilution a) 1:10 + 90 µL Standard 0 (mix thoroughly).

## 6 ASSAY PROCEDURE

### 6.1 General Remarks

- All reagents and specimens must be allowed to come to room temperature before use. All reagents must be mixed without foaming.
- Once the test has been started, all steps should be completed without interruption.
- Use new disposal plastic pipette tips for each standard, control or sample in order to avoid cross contamination.
- Optical density is a function of the incubation time and temperature. Before starting the assay, it is recommended that all reagents are ready, caps removed, all needed wells secured in holder, etc. This will ensure equal elapsed time for each pipetting step without interruption.
- As a general rule the enzymatic reaction is linearly proportional to time and temperature.

## 6.2 Test Procedure

Each run must include a standard curve.

1. Secure the desired number of Microtiter wells in the frame holder.
2. Dispense **25 µL** of each **Standard, Control** and **sample** with new disposable tips into appropriate wells.
3. Dispense **100 µL Enzyme Conjugate** into each well.  
Thoroughly mix for 10 seconds. It is important to have a complete mixing in this step.
4. Incubate for **60 minutes** at room temperature.
5. Rinse the wells **3 times** with **400 µL** diluted Wash Solution per well, if a plate washer is used.  
- OR -  
Briskly shake out the contents of the wells.  
Rinse the wells **3 times** with **300 µL** diluted Wash Solution per well for manual washing. Strike the wells sharply on absorbent paper to remove residual droplets.  
**Important note:**  
The sensitivity and precision of this assay is markedly influenced by the correct performance of the washing procedure!
6. Add **100 µL** of **Substrate Solution** to each well.
7. Incubate for **15 minutes** at room temperature.
8. Stop the enzymatic reaction by adding **50 µL** of **Stop Solution** to each well.
9. Determine the optical density of the solution in each well at **450 nm (reading) and at 620 nm to 630 nm (background subtraction, recommended)** with a microtiter plate reader.  
It is recommended that the wells be read **within 10 minutes** after adding the Stop Solution.

## 6.3 Calculation of Results

1. Calculate the average optical density (OD) values for each set of standards, controls and patient samples.
2. Using graph paper, construct a standard curve by plotting the mean OD obtained from each standard against its concentration with OD value on the vertical (Y) axis and concentration on the horizontal (X) axis.
3. Using the mean OD value for each sample determine the corresponding concentration from the standard curve.
4. Automated method: The results in the Instructions for Use have been calculated automatically using a 4-Parameter curve fit. (4-Parameter Rodbard or 4-Parameter Marquardt are the preferred methods.) Other data reduction functions may give slightly different results.
5. The concentration of the samples can be read directly from this standard curve. Samples with concentrations higher than that of the highest standard have to be further diluted or reported as > 20 ng/mL. For the calculation of the concentrations this dilution factor has to be taken into account.

### 6.3.1 Example of Typical Standard Curve

The following data is for demonstration only and **cannot** be used in place of data generations at the time of assay.

Standard	Optical Density (450 nm)
Standard 0 (0.0 ng/mL)	0.03
Standard 1 (0.4 ng/mL)	0.08
Standard 2 (1.0 ng/mL)	0.16
Standard 3 (4.0 ng/mL)	0.55
Standard 4 (10 ng/mL)	1.27
Standard 5 (20 ng/mL)	2.28

## 7 EXPECTED NORMAL VALUES

It is strongly recommended that each laboratory should determine its own normal and abnormal values. In a study conducted with apparently healthy adults, using the DEMEDITEC AMH ELISA the following data were observed:

Population	n	Mean (ng/mL)	Median (ng/mL)	2.5 <sup>th</sup> - 97.5 <sup>th</sup> Percentile (ng/mL)	Range (min. - max.) (ng/mL)
<b>Males</b>	30	3.99	3.69	0.32 - 7.36	0.06 - 7.75
<b>Females (20 - 29 years)</b>	30	2.71	2.38	0.69 - 6.23	0.66 - 6.37
<b>Females (30 - 39 years)</b>	30	2.42	1.85	0.51 - 6.72	0.48 - 8.39
<b>Females (40 - 49 years)</b>	30	0.61	0.29	< 0.06 - 4.09	< 0.06 - 4.37

Median AMH values in women diagnosed with polycystic ovary syndrome (PCOS) were two- to four-fold higher compared to healthy women and did not markedly decline with increasing age. In women diagnosed with PCOS, a median of 7.19 ng/ml was observed in subjects aged 20–24 years, and a median of 6.46 ng/ml was observed in subjects aged 35–39 years. (11) For women with AMH  $\leq$  0.681 ng/mL, the probability of a low antral follicle count (AFC 0-7) is 63%, the probability of belonging to the middle AFC group (8-15) is approximately 32%, and the probability of AFC > 15 is only 4.4%. (12, 13) The results alone should not be the only reason for any therapeutic consequences. The results should be correlated to other clinical observations and diagnostic tests.

## 8 QUALITY CONTROL

Good laboratory practice requires that controls be run with each calibration curve. A statistically significant number of controls should be assayed to establish mean values and acceptable ranges to assure proper performance. It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results. Use controls at both normal and pathological levels. The controls and the corresponding results of the QC-Laboratory are stated in the QC certificate added to the kit. The values and ranges stated on the QC sheet always refer to the current kit lot and should be used for direct comparison of the results.

It is also recommended to make use of national or international Quality Assessment programs in order to ensure the accuracy of the results. Employ appropriate statistical methods for analysing control values and trends. If the results of the assay do not fit to the established acceptable ranges of control materials patient results should be considered invalid. In this case, please check the following technical areas: Pipetting and timing devices; photometer, expiration dates of reagents, storage and incubation conditions, aspiration and washing methods. After checking the above mentioned items without finding any error contact your distributor or DEMEDITEC directly.

## 9 PERFORMANCE CHARACTERISTICS

### 9.1 Assay Dynamic Range

The range of the assay is between 0.062 ng/mL – 20.0 ng/mL.

### 9.2 Specificity of Antibodies (Cross-Reactivity)

The following substances were tested for cross-reactivity of the assay:

Substance	Added conc. (ng/mL)	Mean cross-reactivity (%)
AMH	0.40 – 10	100
Inhibin A	2.0 – 2000	0.03
Actvin AB	2.0 – 2000	0.27
LH	2.0 – 2000	0.18
FSH	2.0 – 2000	0.26
HCG	2.0 – 2000	0.12
TSH	2.0 – 2000	0.39
TGF- $\beta$ 1	2.0 – 2000	0.18
TGF- $\beta$ 2	2.0 – 2000	0.18
Prolactin	2.0 – 2000	0.38

No substantial cross-reactivity of the assay to structurally related substances is detected.

### 9.3 Sensitivity

The analytical sensitivity of the DEMEDITEC ELISA was calculated by adding 2 standard deviations to the mean of 20 replicate analyses of the Standard 0 and was found to be 0.044 ng/mL.

The Limit of Blank (LoB) is 0.044 ng/mL.

The Limit of Detection (LoD) is 0.052 ng/mL.

The Limit of Quantification (LoQ) is 0.062 ng/mL.

### 9.4 Reproducibility

#### 9.4.1 Intra-Assay

The within-assay variability was determined by measuring each sample 10 times per run (n = 10):

Sample	n	Mean (ng/mL)	CV (%)
1	10	0.22	7.3
2	10	0.67	3.1
3	10	5.76	2.4
4	10	15.88	2.6

#### 9.4.2 Inter-Assay

The between-assay variability was determined by measuring each sample 10 times per run for 3 days (n = 30):

Sample	n	Mean (ng/mL)	CV (%)
1	30	0.23	6.0
2	30	0.69	5.4
3	30	5.73	3.1
4	30	15.90	4.1

### 9.5 Recovery

Samples have been spiked by adding AMH solutions with known concentrations.

The recovery (%) was calculated by multiplying the ratio of measured and expected values with 100.

		Sample 1	Sample 2	Sample 3	Sample 4
Concentration (ng/mL)		0.20	0.70	5.85	15.89
Average Recovery (%)		98.8	97.5	97.0	98.5
Range of Recovery (%)	from	94.2	95.4	94.1	95.9
	to	101.8	100.3	101.1	100.8

### 9.6 Linearity

Samples were measured undiluted and in serial dilutions with standard 0. The recovery (%) was calculated by multiplying the ratio of expected and measured values with 100.


		Sample 1	Sample 2	Sample 3	Sample 4
Concentration (ng/mL)		3.01	5.99	9.75	16.29
Average Recovery (%)		105.2	106.8	95.2	110.3
Range of Recovery (%)	from	94.7	102.4	94.5	106.6
	to	111.9	110.9	95.9	113.8

### 9.7 Comparison Studies

A comparison of the DEMEDITEC AMH ELISA (DE6141) (y) and a reference method (x) using clinical samples gave the following correlation:

$$n = 52$$

$$r = 0.980$$
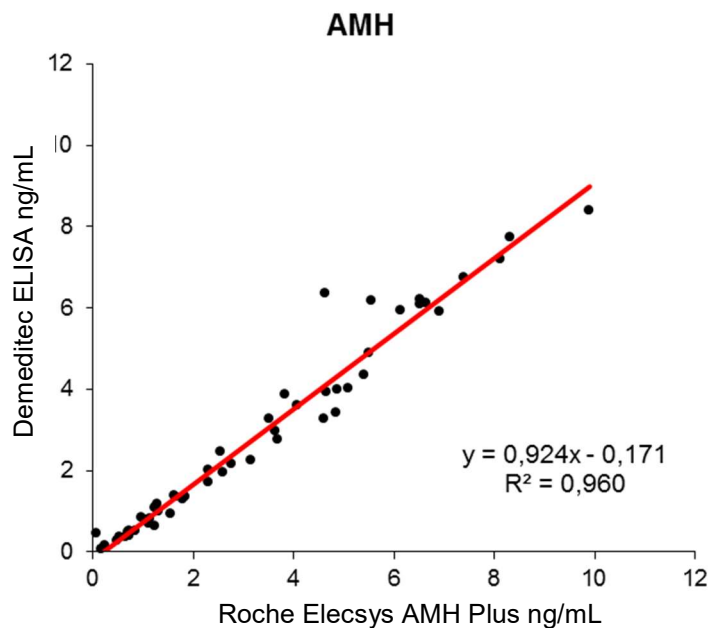
## **10 LIMITATIONS OF USE**

Reliable and reproducible results will be obtained when the assay procedure is performed with a complete understanding of the package insert instruction and with adherence to good laboratory practice.

Any improper handling of samples or modification of this test might influence the results.

### **10.1 Interfering Substances**

Haemoglobin (up to 4 mg/mL), Bilirubin (up to 0.5 mg/mL) and Triglyceride (up to 7.5 mg/mL) have no influence on the assay results.

### **10.2 Drug Interferences**

Until today no substances (drugs) are known to us, which have a substantial influence to the measurement of AMH in a sample.

### **10.3 High-Dose-Hook Effect**

Hook effect was not observed in this test up to a concentration of 400 ng/mL of AMH.

## **11 LEGAL ASPECTS**

### **11.1 Reliability of Results**

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable national standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test. The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications. In case of any doubt or concern please contact DEMEDITEC.

### **11.2 Therapeutic Consequences**

Therapeutic consequences should never be based on laboratory results alone even if all test results are in agreement with the items as stated under point 11.1. Any laboratory result is only a part of the total clinical picture of a patient. Only in cases where the laboratory results are in acceptable agreement with the overall clinical picture of the patient should therapeutic consequences be derived. The test result itself should never be the sole determinant for deriving any therapeutic consequences.

### **11.3 Liability**

Any modification of the test kit and/or exchange or mixture of any components of different lots from one test kit to another could negatively affect the intended results and validity of the overall test. Such modification and/or exchanges invalidate any claim for replacement. Claims submitted due to customer misinterpretation of laboratory results subject to point 11.2 are also invalid. Regardless, in the event of any claim, the manufacturer's liability is not to exceed the value of the test kit. Any damage caused to the test kit during transportation is not subject to the liability of the manufacturer.

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**12 REFERENCES / LITERATURE**

1. Jopling H, Yates A, Burgoyne N, Hayden K, Chaloner C, Tetlow L. Paediatric Anti-Müllerian Hormone measurement: Male and female reference intervals established using the automated Beckman Coulter Access AMH assay. *Endocrinol Diabetes Metab.* 2018;1(4):e00021.
2. Hampl R, Šnajderová, M, Mardešić T. Antimüllerian Hormone (AMH) Not Only a Marker for Prediction of Ovarian Reserve. *Physiol. Res.* 60: 217-223, 2011.
3. Matuszczak E, Hermanowicz A, Komarowska M, Debek W. Serum AMH in Physiology and Pathology of Male Gonads. *Int J Endocrinol.* 2013:128907.
4. Victoria M, Labrosse J, Krief F, Cédric-Durnerin I, Comtet M, Grynberg M. Anti Müllerian Hormone: More than a biomarker of female reproductive function. *J Gynecol Obstet Hum Reprod.* 2019 Jan;48(1):19-24.
5. Broer S, Broekmans F, Laven J, Fauser B. Anti-Müllerian hormone: ovarian reserve testing and its potential clinical implications. *Hum Reprod Update.* 2014 Sep-Oct;20(5):688-701.
6. Broer S, Eijkemans M, Scheffer G, van Rooij I, de Vet A, Themmen A, Laven J, de Jong F, te Velde E, Fauser B, Broekmans F. Anti-Müllerian Hormone Predicts Menopause: A Long-Term Follow-Up Study in Normoovulatory Women. *J Clin Endocrinol Metab.* 2011 Aug;96(8):2532-9.
7. Visser J, de Jong F, Laven, J, Themmen A. Anti-Müllerian hormone: a new marker for ovarian function. *Reproduction.* 2006 Jan;131(1):1-9.
8. Broer S, van Disseldorp J, Broeze K, Dolleman M, Opmeer B, Bossuyt P, Eijkemans M, Mol B, Broekmans F; IMPORT study group. Added value of ovarian reserve testing on patient characteristics in the prediction of ovarian response and ongoing pregnancy: an individual patient data approach. *Hum Reprod Update.* 2013 Jan-Feb;19(1):26-36.
9. Ebner T, Sommergruber M, Moser M, Shebl O, Schreier-Lechner E, Tews G. Basal level of anti-Müllerian hormone is associated with oocyte quality in stimulated cycles. *Hum Reprod.* 2006 Aug;21(8):2022-6.
10. Sun H, Mao H, Cai J, Zhao Y. Research progress on anti-müllerian hormone clinical applications and immunoassay development. *Frontiers in Laboratory Medicine.* 2018 2(1):14-18.
11. Anckaert E, Öktem M, Thies A, Cohen-Bacrie M, Daan NM, Schiettecatte J, Müller C, Topcu D, Gröning A, Ternaux F, Engel C, Engelmann S, Milczynski C. Multicenter analytical performance evaluation of a fully automated anti-Müllerian hormone assay and reference interval determination. *Clinical biochemistry* 2016, 49(3), 260–267.
12. Anderson R A, Anckaert E, Bosch E, Dewailly D, Dunlop CE, Fehr D, Nardo L, Smitz J, Tremellen K, Denk B, Geistanger A, Hund M. Prospective study into the value of the automated Elecsys antimüllerian hormone assay for the assessment of the ovarian growing follicle pool. *Fertility and sterility* 2015, 103(4), 1074–1080.e4.
13. Practice Committee of the American Society for Reproductive Medicine. Testing and interpreting measures of ovarian reserve: a committee opinion. *Fertility and sterility* 2012, 98(6), 1407–1415.

**SYMBOLS USED WITH DEMEDITEC ASSAYS**

Symbol	English	Deutsch	Français	Espanol	Italiano
	European Conformity	CE-Konformitätskennzeichnung	Conforme aux normes européennes	Conformidad europea	Conformità europea
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las Instrucciones	Consultare le istruzioni per l'uso
	In vitro diagnostic device	In-vitro-Diagnostikum	Usage Diagnostic in vitro	Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Référence	Número de catálogo	No. di Cat.
	Lot. No. / Batch code	Chargen-Nr.	No. de lot	Número de lote	Lotto no
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Note warnings and precautions	Warnhinweise und Vorsichtsmaßnahmen beachten	Avertissements et mesures de précaution font attention	Tiene en cuenta advertencias y precauciones	Annoti avvisi e le precauzioni
	Storage Temperature	Lagerungstemperatur	Temperature de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeitsdatum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
<i>Distributed by</i>	Distributor	Vertreiber	Distributeur	Distribuidor	Distributore