
# Product information

Information about other products is available at: [www.demeditec.com](http://www.demeditec.com)

## CEA ELISA

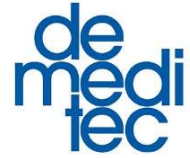**DE1868**

**96**

Demeditec Diagnostics GmbH  
Lise-Meitner-Strasse 2  
24145 Kiel – Germany  
[www.demeditec.com](http://www.demeditec.com)

**Please use only the valid version of the Instructions for Use provided with the kit.  
Verwenden Sie nur die jeweils gültige, im Testkit enthaltene, Gebrauchsanweisung.  
Si prega di usare la versione valida delle istruzioni per l'uso a disposizione con il kit.  
Por favor, se usa solo la version valida de la metodico técnico incluido aqui en el kit.  
Utilisez seulement la version valide des Instructions d'utilisation fournies avec le kit.**

**Table of Contents / Inhaltsverzeichnis / Tabella die Contenuti / Tabla de Contenidos / Sommaire**

1	INTRODUCTION.....	4
2	PRINCIPLE OF THE TEST.....	4
3	WARNINGS AND PRECAUTIONS.....	5
4	REAGENTS.....	6
5	SPECIMEN COLLECTION AND PREPARATION.....	7
6	ASSAY PROCEDURE.....	7
7	EXPECTED NORMAL VALUES.....	8
8	QUALITY CONTROL.....	9
9	PERFORMANCE CHARACTERISTICS.....	9
10	LIMITATIONS OF USE.....	10
11	LEGAL ASPECTS.....	10
1	EINLEITUNG.....	11
2	TESTPRINZIP.....	11
3	VORSICHTSMAßNAHMEN.....	11
4	BESTANDTEILE DES KITS.....	12
5	PROBENVORBEREITUNG.....	13
6	TESTDURCHFÜHRUNG.....	13
7	ERWARTETE WERTE.....	15
8	QUALITÄTSKONTROLLE.....	15
9	ASSAY-CHARAKTERISTIKA.....	15
10	GRENZEN DES TESTS.....	15
11	RECHTLICHE GRUNDLAGEN.....	16

---

1	INTRODUZIONE .....	17
2	PRINCIPIO DEL TEST .....	17
3	PRECAUZIONI .....	17
4	COMPONENTI DEL KIT .....	18
5	CAMPIONI.....	19
6	ATTUAZIONE DEL TEST .....	19
7	VALORI NORMALI.....	21
8	CONTROLLO QUALITÀ.....	21
9	CARATTERISTICHE DEL TEST.....	21
10	LIMITAZIONE DEL TEST.....	22
11	ASPETTI LEGALI .....	22
1	INTRODUCCIÓN .....	23
2	FUNDAMENTO DEL ENSAYO .....	23
3	PRECAUCIONES.....	23
4	COMPONENTES DEL KIT.....	24
5	MUESTRAS.....	25
6	PROCEDIMIENTO DE ENSAYO .....	25
7	VALORES ESPERADOS .....	27
8	CONTROL DE CALIDAD .....	27
9	CARACTERÍSTICAS DEL ENSAYO .....	27
10	LIMITACIONES DE USO .....	28
11	ASPECTOS LEGALES .....	28
1	INTRODUCTION.....	29
2	PRINCIPE DU TEST .....	29
3	PRECAUTIONS D'UTILISATION.....	29
4	COMPOSITION DU KIT .....	30
5	ECHANTILLON .....	31
6	RÉALISATION DU TEST .....	32
7	VALEURS ATTENDUES.....	33
8	CONTROLE DE QUALITE .....	33
9	CARACTERISTIQUES DU TEST .....	33
10	LIMITES D'UTILISATION.....	34
11	ASPECTS LEGAUX .....	34
12	REFERENCES / LITERATURE .....	35
	SYMBOLS USED WITH DEMEDITEC ASSAYS .....	36

## 1 INTRODUCTION

### 1.1 Intended Use

The **Demeditec CEA ELISA** is an enzyme immunoassay for the quantitative *in vitro diagnostic* measurement of Carcinoembryonic Antigen (CEA) in serum or lithium heparin plasma.

### 1.2 Summary and Explanation

Carcinoembryonic antigen (CEA) is a 200 kD glycoprotein. Cells that express CEA, both incorporate this glycoprotein into their cell membrane and release it into the blood. CEA is therefore detectable both on cells and in body liquids.

The normal values are normally < 5 ng/mL

Tumors that are associated with increased CEA values are:

Colon carcinoma, stomach carcinoma, breast cancer, lung cancer, pancreatic cancer, gullet cancer.

The most important role of CEA is in colon cancer, since the level of CEA correlates with the stage of the tumor.

## 2 PRINCIPLE OF THE TEST

The Demeditec CEA ELISA Kit is a solid phase enzyme-linked immunosorbent assay (ELISA) based on the **sandwich principle**. The microtiter wells are coated with a monoclonal antibody directed towards a unique antigenic site of the CEA molecule. An aliquot of patient sample containing endogenous CEA is incubated in the coated well with enzyme conjugate, which is an anti-CEA antibody conjugated with horseradish peroxidase. After incubation the unbound conjugate is washed off.

The amount of bound peroxidase conjugate is proportional to the concentration of CEA in the sample. Having added the substrate solution, the intensity of colour developed is proportional to the concentration of CEA in the patient sample.

### 3 WARNINGS AND PRECAUTIONS

- This kit is for in vitro diagnostic use only. For professional use only.
- All reagents of this test kit which contain human serum or plasma have been tested and confirmed negative for HIV I/II, HBsAg and HCV by FDA approved procedures. All reagents, however, should be treated as potential biohazards in use and for disposal.
- Before starting the assay, read the instructions completely and carefully. Use the valid version of instructions for use provided with the kit. Be sure that everything is understood.
- The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch and used in the frame provided.
- Pipetting of samples and reagents must be done as quickly as possible and in the same sequence for each step.
- Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution colored. Do not pour reagents back into vials as reagent contamination may occur.
- Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells.
- Do not let wells dry during assay; add reagents immediately after completing the rinsing steps.
- Allow the reagents to reach room temperature (21 °C to 26 °C) before starting the test. Temperature will affect the absorbance readings of the assay. However, values for the patient samples will not be affected.
- Never pipet by mouth and avoid contact of reagents and specimens with skin and mucous membranes.
- Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.
- Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents or specimens may give false results.
- Handling should be done in accordance with the procedures defined by an appropriate national biohazard safety guideline or regulation.
- Do not use reagents beyond expiry date as shown on the kit labels.
- All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microtiterplate readers.
- Do not mix or use components from kits with different lot numbers. It is advised not to exchange wells of different plates even of the same lot. The kits may have been shipped or stored under different conditions and the binding characteristics of the plates may result slightly different.
- Avoid contact with *Stop Solution* containing 0.5 M H<sub>2</sub>SO<sub>4</sub>. It may cause skin irritation and burns.
- Some reagents contain Proclin 300, BND and/or MIT as preservatives. In case of contact with eyes or skin, flush immediately with water.
- TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them. If inhaled, take the person to open air.
- Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guideline or regulation.
- For information on hazardous substances included in the kit please refer to Safety Data Sheets. Safety Data Sheets for this product are available upon request directly from Demeditec.

## 4 REAGENTS

### 4.1 Reagents provided

1. **SORB MT Microtiterwells**, 12 x 8 (break apart) strips, 96 wells; Wells coated with anti-CEA antibody (monoclonal).
2. **CAL 0 Zero Standard**, 1 vial, 3 mL, ready to use; Contains non-mercury preservative.
3. **CAL 1 – 5 Standard (Standard 1-5)**, 5 vials, 1 mL, ready to use; Concentrations: 5 - 10 - 25 - 50 - 100 ng/mL; Contain non-mercury preservative.
4. **CONTROL low & high LYO Control Low & High**, 2 vials, (lyophilized), 1 mL each, See “Reagent Preparation”. For control values and ranges please refer to QC-Datasheet. Contain non-mercury preservative.
5. **ENZ CONJ Enzyme Conjugate**, 1 vial, 14 mL, ready to use, Monoclonal anti-CEA antibody conjugated to horseradish peroxidase; Contains non-mercury preservative.
6. **SUB TMB Substrate Solution**, 1 vial, 14 mL, ready to use, Tetramethylbenzidine (TMB).
7. **STOP SOLN Stop Solution**, 1 vial, 14 mL, ready to use, contains 0.5 M H<sub>2</sub>SO<sub>4</sub>. Avoid contact with the stop solution. It may cause skin irritations and burns.
8. **WASH SOLN 40x Wash Solution**, 1 vial, 30 mL (40X concentrated), See “Reagent Preparation”.

**Note:** Additional *Zero Standard* for sample dilution is available upon request.

### 4.2 Materials required but not provided

- A microtiter plate calibrated reader (450 nm, with reference wavelength at 620 nm - 630 nm)
- Calibrated variable precision micropipettes
- Absorbent paper
- Distilled or deionized water
- Timer
- Linear graph paper or software for data reduction

### 4.3 Storage Conditions

When stored at 2-8 °C unopened reagents will retain reactivity until expiration date. Do not use reagents beyond this date. Opened reagents must be stored at 2-8 °C. Microtiter wells must be stored at 2-8 °C. Once the foil bag has been opened, care should be taken to close it tightly again. Opened kits retain activity for two months if stored as described above.

### 4.4 Reagent Preparation

Bring all reagents and required number of strips to room temperature prior to use.

#### Control

Reconstitute the lyophilized content with 1.0 mL deionized water and let stand for 10 minutes in minimum. Mix the control several times before use.

**Note:** *The reconstituted control should be apportioned and stored at -20 °C.*

#### Wash Solution

Add deionized water to the 40X concentrated Wash Solution.

Dilute 30 mL of concentrated *Wash Solution* with 1170 mL deionized water to a final volume of 1200 mL. *The diluted Wash Solution is stable for 2 weeks at room temperature.*

### 4.5 Disposal of the Kit

The disposal of the kit must be made according to the national regulations. Special information for this product is given in the Safety Data Sheet.

### 4.6 Damaged Test Kits

In case of any severe damage to the test kit or components, Demeditec has to be informed in writing, at the latest, one week after receiving the kit. Severely damaged single components should not be used for a test run. They have to be stored until a final solution has been found. After this, they should be disposed according to the official regulations.

## 5 SPECIMEN COLLECTION AND PREPARATION

Serum or lithium heparin plasma can be used in this assay.

Do not use haemolytic, icteric or lipaemic specimens.

*Note:* Samples containing sodium azide should not be used in the assay.

### 5.1 Specimen Collection

#### Serum:

Collect blood by venipuncture (e.g. Sarstedt Monovette for serum), allow to clot, and separate serum by centrifugation at room temperature. Do not centrifuge before complete clotting has occurred. Patients receiving anticoagulant therapy may require increased clotting time.

#### Plasma:

Whole blood should be collected into centrifuge tubes containing anti-coagulant (e.g. Sarstedt Monovette with the appropriate plasma preparation) and centrifuged immediately after collection.

### 5.2 Specimen Storage and Preparation

Specimens should be capped and may be stored for up to 48 hours at 2-8 °C prior to assaying.

Specimens held for a longer time (up to 12 months) should be frozen only once at -20 °C prior to assay. Thawed samples should be inverted several times prior to testing.

### 5.3 Specimen Dilution

If in an initial assay, a specimen is found to contain more than the highest standard, the specimens can be diluted with *Zero Standard* and reassayed as described in Assay Procedure.

For the calculation of the concentrations this dilution factor has to be taken into account.

#### Example:

- a) dilution 1:10:            10 µL sample + 90 µL *Zero Standard* (mix thoroughly)  
b) dilution 1:100:        10 µL dilution a) 1:10 + 90 µL *Zero Standard* (mix thoroughly).

## 6 ASSAY PROCEDURE

### 6.1 General Remarks

- All reagents and specimens must be allowed to come to room temperature before use. All reagents must be mixed without foaming.
- Once the test has been started, all steps should be completed without interruption.
- Use new disposal plastic pipette tips for each standard, control or sample in order to avoid cross contamination.
- Absorbance is a function of the incubation time and temperature. Before starting the assay, it is recommended that all reagents are ready, caps removed, all needed wells secured in holder, etc. This will ensure equal elapsed time for each pipetting step without interruption.
- As a general rule the enzymatic reaction is linearly proportional to time and temperature.

## 6.2 Test Procedure

Each run must include a standard curve.

All standards, samples, and controls should be run in duplicate. All standards, samples, and controls should be run concurrently so that all conditions of testing are the same.

1. Secure the desired number of Microtiter wells in the frame holder.
2. Dispense **50 µL** of each **Standard, Control** and **samples** with new disposable tips into appropriate wells.
3. Dispense **100 µL Enzyme Conjugate** into each well.  
Thoroughly mix for 10 seconds. It is important to have a complete mixing in this step.
4. Incubate for **60 minutes** at room temperature.
5. Briskly shake out the contents of the wells.  
Rinse the wells **3 times** with 400 µL diluted *Wash Solution* per well, if a plate washer is used - or - rinse the wells 3 times with 300 µL diluted *Wash Solution* per well for manual washing.  
Strike the wells sharply on absorbent paper to remove residual droplets.  
**Important note:** The sensitivity and precision of this assay is markedly influenced by the correct performance of the washing procedure!
6. Add **100 µL** of **Substrate Solution** to each well.
7. Incubate for **30 minutes** at room temperature.
8. Stop the enzymatic reaction by adding **100 µL** of **Stop Solution** to each well.
9. Determine the absorbance (OD) of the solution in each well at 450 nm (reading) and at 620 - 630 nm (background subtraction, recommended) with a microtiter plate reader.  
It is recommended that the wells be read **within 10 minutes** after adding the *Stop Solution*.

## 6.3 Calculation of Results

1. Calculate the average absorbance values for each set of standards, controls and patient samples.
2. Using linear graph paper, construct a standard curve by plotting the mean absorbance obtained from each standard against its concentration with absorbance value on the vertical (Y) axis and concentration on the horizontal (X) axis.
3. Using the mean absorbance value for each sample determine the corresponding concentration from the standard curve.
4. Automated method: The results in the Instructions for Use have been calculated automatically using a 4-Parameter curve fit. (4 Parameter Rodbard or 4 Parameter Marquardt are the preferred methods.)  
Other data reduction functions may give slightly different results.
5. The concentration of the samples can be read directly from this standard curve. Samples with concentrations higher than that of the highest standard have to be further diluted or reported as > 100 ng/mL. For the calculation of the concentrations this dilution factor has to be taken into account.

### 6.3.1 Example of Typical Standard Curve

The following data is for demonstration only and **cannot** be used in place of data generations at the time of assay.

Standard	Optical Units (450 nm)
Zero Standard (0 ng/mL)	0.06
Standard 1 (5 ng/mL)	0.20
Standard 2 (10 ng/mL)	0.34
Standard 3 (25 ng/mL)	0.62
Standard 4 (50 ng/mL)	1.12
Standard 5 (100 ng/mL)	2.04

## 7 EXPECTED NORMAL VALUES

It is strongly recommended that each laboratory should determine its own normal and abnormal values. The literature cut-off is 5 ng/mL for non-smokers and 10 ng/mL for smokers. The results alone should not be the only reason for any therapeutic consequences. The results should be correlated to other clinical observations and diagnostic tests.

## 8 QUALITY CONTROL

Good laboratory practice requires that controls be run with each calibration curve. A statistically significant number of controls should be assayed to establish mean values and acceptable ranges to assure proper performance. It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results. Use controls at both normal and pathological levels. The controls and the corresponding results of the QC-Laboratory are stated in the QC certificate added to the kit. The values and ranges stated on the QC sheet always refer to the current kit lot and should be used for direct comparison of the results.

It is also recommended to make use of national or international Quality Assessment programs in order to ensure the accuracy of the results. Employ appropriate statistical methods for analysing control values and trends. If the results of the assay do not fit to the established acceptable ranges of control materials patient results should be considered invalid. In this case, please check the following technical areas: Pipetting and timing devices; photometer, expiration dates of reagents, storage and incubation conditions, aspiration and washing methods. After checking the above mentioned items without finding any error contact your distributor or Demeditec directly.

## 9 PERFORMANCE CHARACTERISTICS

### 9.1 Assay Dynamic Range

The range of the assay is between 0.596 ng/mL – 100 ng/mL.

### 9.2 Specificity of Antibodies (Cross Reactivity)

Cross-reactivities are not known.

### 9.3 Sensitivity

The analytical sensitivity of the Demeditec ELISA was calculated by adding 2 standard deviations to the mean of 20 replicate analyses of the *Zero Standard* and was found to be < 0.596 ng/mL.

### 9.4 Reproducibility

Intra-assay variations of 3.2% – 4.9% were observed.

Inter-assay variations of 5.6% – 6.3% were observed.

### 9.5 Recovery

Samples have been spiked by adding CEA solutions with known concentrations in a 1:1 ratio.

The % recovery has been calculated by multiplication of the ratio of the measurements and the expected values with 100 (expected value = (endogenous CEA + added CEA) / 2; because of a 1:2 dilution of serum with spike material).

	Sample 1	Sample 2	Sample 3
<b>Concentration [ng/mL]</b>	3.3	18.3	47.1
<b>Average Recovery [%]</b>	90.3	89.4	101.1
<b>Range of Recovery [%]</b>	from	87.6	85.8
	to	95.6	95.2
			107.0

### 9.6 Linearity

	Sample 1	Sample 2	Sample 3
<b>Concentration [ng/mL]</b>	17.7	44.6	85.3
<b>Average Recovery [%]</b>	109.3	96.7	99.9
<b>Range of Recovery [%]</b>	from	105.7	89.9
	to	114.6	100.4
			106.3

## **10 LIMITATIONS OF USE**

Reliable and reproducible results will be obtained when the assay procedure is performed with a complete understanding of the package insert instruction and with adherence to good laboratory practice. Any improper handling of samples or modification of this test might influence the results.

### **10.1 Interfering Substances**

Haemolytic, icteric and lipemic sera should be avoided. The assay contains reagents to minimize interference of HAMA and heterophilic antibodies. However, extremely high titers of HAMA or heterophilic antibodies may interfere with the test results.

### **10.2 Drug Interferences**

Smokers have been reported to exhibit increased levels of CEA (see expected values)

### **10.3 High-Dose-Hook Effect**

No hook effect was observed in this test up to 8000 ng/mL of CEA.

## **11 LEGAL ASPECTS**

### **11.1 Reliability of Results**

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable national standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test. The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications. In case of any doubt or concern please contact Demeditec.

### **11.2 Therapeutic Consequences**

Therapeutic consequences should never be based on laboratory results alone even if all test results are in agreement with the items as stated under point 11.1. Any laboratory result is only a part of the total clinical picture of a patient. Only in cases where the laboratory results are in acceptable agreement with the overall clinical picture of the patient should therapeutic consequences be derived. The test result itself should never be the sole determinant for deriving any therapeutic consequences.

### **11.3 Liability**

Any modification of the test kit and/or exchange or mixture of any components of different lots from one test kit to another could negatively affect the intended results and validity of the overall test. Such modification and/or exchanges invalidate any claim for replacement. Claims submitted due to customer misinterpretation of laboratory results subject to point 11.2 are also invalid. Regardless, in the event of any claim, the manufacturer's liability is not to exceed the value of the test kit. Any damage caused to the test kit during transportation is not subject to the liability of the manufacturer.

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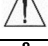
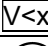
## 12 REFERENCES / LITERATURE

1. Engvall, E. Methods in Enzymology, Volume 70, Van Vunakis, H. and Langone, J.J. (eds.), Academic Press, New York, 419-492 (1980).
2. Uotila, M. Ruoslahti, E. and Engvall, E., J. Immunol. Methods, 42, 11-15 (1981).
3. Gold, P., Freedman S.O: Demonstration of tumor specific antigen in human colonic carcinoma by immunologic tolerance and absorption techniques. J Exp. Med. 1965; 127:439-462.
4. Thompson, D.P.M. Krupey J, Freedman S.O., et al. The radioimmunoassay of circulating carcinoembryonic antigen of the human digestive system. Proc Natl Acad Sci USA 1969, 64:161-167.
5. Schwartz M.K. Tumor markers in diagnosis and screening. In: Ting S.W., Chen J.S., Schwartz M.K: , eds. Human tumor markers, Amsterdam: Elsevier Science, 1987; 3-16.
6. Zamchek N, and Martin E.W. Sequential Carcinoembryonic Antigen Levels in Pancreatic Cancer: Some Clinical Correlations. Cancer 1981; 47:1620-1627.
7. Mughal A.W., Hortobagyi G.N.; Fritsche H.A., Buzdar A.U. Yap H-Y. and Blumenschein G.R. Serial Plasma Carcinoembryonic Antigen Measurements During Treatment of Metastatic Breast Cancer. JAMA 1983; 259; 1881-1886.

## SYMBOLS USED WITH DEMEDITEC ASSAYS

Symbol	English	Deutsch	Française	Espanol	Italiano
	European Conformity	CE-Konformitätskennzeichnung	Conforme aux normes européennes	Conformidad europea	Conformità europea
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las Instrucciones	Consultare le istruzioni per l'uso
	In vitro diagnostic device	In-vitro-Diagnostikum	utilisation Diagnostic in vitro	Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Référence	Número de catálogo	No. di catalogo
	Lot. No. / Batch code	Chargen-Nr.	No. de lot	Número de lote	Lotto no
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Note warnings and precautions	Warnhinweise und Vorsichtsmaßnahmen beachten	Avertissements et mesures de précaution font attention	Tiene en cuenta y advertencias y precauciones	Annoti avvisi e le precauzioni
	Storage Temperature	Lagerungstemperatur	Température de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeitsdatum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
<i>Distributed by</i>	Distributed by	Vertrieb durch	Distribution par	Distribución por	Distribuzione da parte di
	Version	Version	Version	Versión	Versione
	Single-use	Einmalverwendung	À usage unique	Uso único	Uso una volta