
# Product information

Information about other products is available at: [www.demeditec.com](http://www.demeditec.com)

---
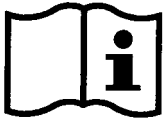
# Parvovirus B19 IgG ELISA

**DE3503**

**96**

Demeditec Diagnostics GmbH  
Lise-Meitner-Strasse 2  
24145 Kiel – Germany  
[www.demeditec.com](http://www.demeditec.com)

**Please use only the valid version of the Instructions for Use provided with the kit.  
Verwenden Sie nur die jeweils gültige, im Testkit enthaltene, Gebrauchsanweisung.  
Si prega di usare la versione valida delle istruzioni per l'uso a disposizione con il kit.  
Por favor, se usa solo la version valida de la metodico técnico incluido aqui en el kit.  
Utilisez seulement la version valide des Instructions d'utilisation fournies avec le kit.**

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## 1 INTENDED USE

The DEMEDITEC Parvovirus B19 IgG (Recombinant) Enzyme Immunoassay Kit provides materials for the qualitative and semi-quantitative determination of IgG-class antibodies to Parvovirus B19 in human serum and plasma (EDTA, lithium heparin or citrate plasma).

**This assay is intended for in vitro diagnostic use only.**

### 1.1 Summary and Explanation

Parvoviruses are cubic single-stranded DNA viruses of about 18-32 nm lacking an envelope. Parvovirus B19 infects only humans, and since there are no cross-reactivities between animal parvoviruses and B19, transmission between pets and humans is not possible. Parvovirus B19 is the causative agent of Erythema infectiosum, the so-called "fifth disease", a mild rash illness that occurs most commonly in children. Infected persons are contagious during the early part of the illness before the rash appears so in adults the rate of epidemia amounts to about 60%. About 20% of adults and children who are infected with parvovirus B19 do not develop any symptoms. Persons infected with the virus, however, do develop lasting immunity that protects them against infection in the future. Parvovirus B19 infection may cause a serious illness in persons with sickle-cell disease or similar types of chronic anemia as well as in persons who have problems with their immune system (people with leukemia or cancer, who are born with immune deficiencies, who have received an organ transplant, or who have HIV infection). Occasionally (less than 5% of all pregnant women infected with parvovirus B19) serious complications may develop during pregnancy: risk of Morbus haemolyticus fetalis.

## 2 PRINCIPLE OF THE TEST

The **Parvovirus B19 IgG ELISA** is a solid phase enzyme-linked immunosorbent assay (ELISA) Microtiter wells as a solid phase are coated with recombinant Parvovirus B19 antigen (VP1-s proteins). **Diluted patient** specimens and **ready-for-use controls** are pipetted into these wells. During incubation Parvovirus B19-specific antibodies of positive specimens and controls are bound to the immobilized antigens. After a washing step to remove unbound sample and control material horseradish peroxidase conjugated anti-human IgG antibodies are dispensed into the wells. During a second incubation this anti-IgG conjugate binds specifically to IgG antibodies resulting in the formation of enzyme-linked immune complexes. After a second washing step to remove unbound conjugate the immune complexes formed (in case of positive results) are detected by incubation with TMB substrate and development of a blue color. The blue color turns into yellow by stopping the enzymatic indicator reaction with sulfuric acid. The intensity of this color is directly proportional to the amount of Parvovirus B19-specific IgG antibody in the patient specimen. Optical density at 450 nm is read using an ELISA microtiter plate reader.

### 3 WARNINGS AND PRECAUTIONS

- This kit is for in vitro diagnostic use only. For professional use only.
- Before starting the assay, read the instructions completely and carefully. Use the valid version of the package insert provided with the kit. Be sure that everything is understood.
- All reagents of this test kit which contain human serum or plasma have been tested and confirmed negative for HIV I/II, HBsAg and HCV by FDA approved procedures. All reagents, however, should be treated as potential biohazards in use and for disposal.
- Avoid contact with Stop Solution containing 0.2 mol/L H<sub>2</sub>SO<sub>4</sub>. It may cause skin irritation and burns.
- TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them. If inhaled, take the person to open air.
- The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch and used in the frame provided
- Pipetting of samples and reagents must be done as quickly as possible and in the same sequence for each step.
- Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution colored. Do not pour reagents back into vials as reagent contamination may occur.
- Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells.
- Do not let wells dry during assay; add reagents immediately after completing the rinsing steps.
- Allow the reagents to reach room temperature (20 °C to 25 °C) before starting the test. Temperature will affect the optical density readings of the assay. However, values for the patient samples will not be affected.
- Never pipette by mouth and avoid contact of reagents and specimens with skin and mucous membranes.
- Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.
- Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents or specimens may give false results.
- Handling should be in accordance with the procedures defined by an appropriate national biohazard safety guideline or regulation.
- Do not use reagents beyond expiry date as shown on the kit labels.
- All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microtiter plate readers.
- Do not mix or use components from kits with different lot numbers. It is advised not to exchange wells of different plates even of the same lot. The kits may have been shipped or stored under different conditions and the binding characteristics of the plates may result slightly different.
- Chemicals and prepared or used reagents have to be treated as hazardous waste according to the national biohazard safety guideline or regulation.
- For information on hazardous substances included in the kit please refer to Safety Data Sheets. Safety Data Sheets for this product are available upon request directly from DEMEDITEC.

## 4 REAGENTS

### 4.1 Reagents provided

1. **SORB MT** Microtiterwells, 12 x 8 (break apart) strips, 96 wells; Wells coated with recombinant Parvovirus B19 antigen (VP1-s protein). (incl. 1 cover foil)
2. **SAM DIL** Sample Diluent \*, 1 vial, 100 mL, ready to use, colored yellow; pH 7.2 ± 0.2.
3. **CAL A** Neg. Control \*, 1 vial, 2.0 mL, ready to use; colored yellow, yellow cap.
4. **CAL B** Cut-off Control \*, 1 vial, 2.0 mL, ready to use; colored yellow, black cap.
5. **CAL C** Pos. Control \*, 1 vial, 1.0 mL, ready to use; colored yellow, red cap.
6. **ENZ CONJ** Enzyme Conjugate \*, 1 vial, 20 mL, ready to use, colored red, antibody to human IgG conjugated to horseradish peroxidase.
7. **SUB TMB** Substrate Solution, 1 vial, 14 mL, ready to use, Tetramethylbenzidine (TMB).
8. **STOP SOLN** Stop Solution, 1 vial, 14 mL, ready to use, contains 0.2 mol/L H<sub>2</sub>SO<sub>4</sub>. Avoid contact with the stop solution. It may cause skin irritations and burns.
9. **WASH SOLN 20x** Wash Solution \*, 1 vial, 30 mL (20X concentrated for 600 mL), pH 6.5 ± 0.1 see „Preparation of Reagents“.

\* contain non-mercury preservative

#### 4.1.1 Material required but not provided

- A calibrated microtiter plate reader (450 nm/620 nm ± 10 nm)
- Calibrated variable precision micropipettes
- Incubator 37 °C
- Manual or automatic equipment for rinsing wells
- Vortex tube mixer
- Freshly distilled water
- Timer
- Absorbent paper

### 4.2 Storage Conditions

When stored at 2 °C to 8 °C unopened reagents will retain reactivity until expiration date. Do not use reagents beyond this date.

Opened reagents must be stored at 2 °C to 8 °C. Microtiter wells must be stored at 2 °C to 8 °C. Once the foil bag has been opened, care should be taken to close it tightly again.

Opened kits retain activity for 8 weeks if stored as described above.

### 4.3 Reagent Preparation

Allow all reagents and required number of strips to reach room temperature prior to use.

#### **Wash Solution**

Dilute *Wash Solution 1+19* (e.g. 10 mL + 190 mL) with fresh and germ free redistilled water. This diluted wash solution has a pH value of 7.2 ± 0.2.

Consumption: ~ 5 mL per determination.

Crystals in the solution disappear by warming up to 37 °C in a water bath. Be sure that the crystals are completely dissolved before use.

*The diluted Wash Solution is stable for 1 week at 2 °C to 8 °C.*

### 4.4 Disposal of the Kit

The disposal of the kit must be made according to the national regulations. Special information for this product is given in the Safety Data Sheets.

### 4.5 Damaged Test Kits

In case of any severe damage to the test kit or components, DEMEDITEC has to be informed in writing, at the latest, one week after receiving the kit. Severely damaged single components should not be used for a test run. They have to be stored until a final solution has been found. After this, they should be disposed according to the official regulations.

## 5 SPECIMEN COLLECTION AND PREPARATION

Serum or plasma (EDTA, lithium heparin or citrate plasma) can be used in this assay.

*Please note:* Samples containing sodium azide should not be used in the assay.

In general, it should be avoided to use haemolytic, icteric or lipaemic specimens. For further information refer to chapter "Interfering Substances".

### 5.1 Specimen Collection

#### Serum:

Collect blood by venipuncture (e.g. Sarstedt Monovette for serum), allow to clot, and separate serum by centrifugation at room temperature. Do not centrifuge before complete clotting has occurred. Patients receiving anticoagulant therapy may require increased clotting time.

#### Plasma:

Whole blood should be collected into centrifuge tubes containing anti-coagulant (e.g. Sarstedt Monovette with the appropriate plasma preparation) and centrifuged immediately after collection.

### 5.2 Specimen Storage

Specimens should be capped and may be stored for up to 5 days at 2 °C to 8 °C prior to assaying.

Specimens held for a longer time (up to 18 months) should be frozen only once at -20 °C prior to assay. Thawed samples should be inverted several times prior to testing.

### 5.3 Specimen Dilution

Prior to assaying dilute each patient specimen **1+100** with *Sample Diluent*;

e.g. 10 µL of specimen + 1 mL of *Sample Diluent*, **mix well, let stand for 15 minutes, mix well again.**

**Please note:** Controls are ready for use and must not be diluted!

## 6 ASSAY PROCEDURE

### 6.1 General Remarks

- **It is very important to bring all reagents, samples and controls to room temperature before starting the test run!**
- Once the test has been started, all steps should be completed without interruption.
- Use new disposal plastic pipette tips for each standard, control or sample in order to avoid cross contamination
- Optical density is a function of the incubation time and temperature. Before starting the assay, it is recommended that all reagents are ready, caps removed, all needed wells secured in holder, etc. This will ensure equal elapsed time for each pipetting step without interruption.
- As a general rule the enzymatic reaction is linearly proportional to time and temperature.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense conjugate without splashing accurately to the bottom of wells.
- During incubation cover microtiter strips with foil to avoid evaporation.

## 6.2 Test Procedure

Prior to commencing the assay, dilute Wash Solution, **prepare patient samples as described in point 5.3**, mix well before pipetting. Carefully establish the **distribution and identification plan** supplied in the kit for all specimens and controls.

1. Select the required number of microtiter strips or wells and insert them into the holder.

Please allocate at least:

- |         |                |                                    |
|---------|----------------|------------------------------------|
| 1 well  | (e.g. A1)      | for the <i>Neg. Control</i> ,      |
| 2 wells | (e.g. B1 + C1) | for the <i>Cut-off Control</i> and |
| 1 well  | (e.g. D1)      | for the <i>Pos. Control</i> .      |

It is left to the user to determine controls and patient samples in duplicate.

2. Dispense  
**100 µL** of *Neg. Control* into well A1  
**100 µL** of *Cut-off Control* into wells B1 + C1  
**100 µL** of *Pos. Control* into well D1 and  
**100 µL** of each diluted sample with new disposable tips into appropriate wells.
3. Cover wells with foil supplied in the kit. Incubate for **60 minutes at 37 °C**.
4. Briskly shake out the contents of the wells.  
Rinse the wells **5 times** with diluted *Wash Solution* (**300 µL per well**). Strike the wells sharply on absorbent paper to remove residual droplets.  
**Important note:** The sensitivity and precision of this assay is markedly influenced by the correct performance of the washing procedure!
5. Dispense **100 µL Enzyme Conjugate** into each well.
6. Incubate for **30 minutes at room temperature (20 °C to 25 °C)**.  
*Do not expose to direct sun light!*
7. Briskly shake out the contents of the wells.  
Rinse the wells **5 times** with diluted *Wash Solution* (300 µL per well). Strike the wells sharply on absorbent paper to remove residual droplets.
8. Add **100 µL** of *Substrate Solution* into all wells.
9. Incubate for **exactly 15 minutes at room temperature (20 °C to 25 °C) in the dark**.  
Any blue color developed during the incubation turns into yellow.  
**Note:** Highly positive patient samples can cause dark precipitates of the chromogen!
11. Read the optical density at **450/620 nm** with a microtiter plate reader **within 30 minutes** after adding the *Stop Solution*.

## 6.3 Measurement

**Measure the optical density (OD)** of all wells **at 450 nm** and record the OD values for each control and patient sample in the distribution and identification plan.

Dual wavelength reading using 620 nm as reference wavelength is recommended.

Where applicable **calculate the mean OD values** of all duplicates.

## 7 RESULTS

### 7.1 Validation of the Test Run

The test run may be considered valid provided the following criteria are met:

<b>Neg. Control</b> in A1:	OD value <b>lower than 0.200</b>
<b>Cut-off Control</b> in B1/C1:	OD value <b>between 0.350 – 0.850</b>
<b>Pos. Control</b> in D1:	OD value <b>between 0.650 – 3.000</b>

**The OD value of the Pos. Control should be higher than the OD value of the Cut-off Control.**

(OD Pos. Control > OD Cut-off Control).

### 7.2 Calculation

#### Mean OD value of Cut-off Control [CO]

Calculate the mean OD value of the two (2) Cut-off Control determinations (e.g. in B1/C1).

**Example:**  $(0.40 + 0.45) / 2 = 0.425 = CO$

### 7.3 Interpretation

<b>POSITIVE</b>	Patient (mean) OD values more than 20% above CO (Mean OD <sub>patient</sub> > 1.2 × CO)
<b>GREY ZONE</b>	Patient (mean) OD values from 20% above to 10% below CO repeat test 2 - 4 weeks later - with <u>new</u> patient samples ( $0.90 \times CO \leq \text{Mean OD}_{\text{patient}} \leq 1.2 \times CO$ )
<b>NEGATIVE</b>	Results in the second test again in the grey zone ⇒ <b>NEGATIVE</b> Patient (mean) OD values more than 10% below CO (Mean OD <sub>patient</sub> < 0.90 × CO)

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## 8 QUALITY CONTROL

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results. Use controls at both normal and pathological levels. It is also recommended to make use of national or international Quality Assessment programs in order to ensure the accuracy of the results. If the results of the assay do not fit to the established acceptable ranges of control materials patient results should be considered invalid. In this case, please check the following technical areas: Pipetting and timing devices; photometer, expiration dates of reagents, storage and incubation conditions, aspiration and washing methods. After checking the above mentioned items without finding any error contact your distributor or DEMEDITEC directly.

## 9 ASSAY CHARACTERISTICS

### 9.1 Analytical Sensitivity

The analytical sensitivity of the DEMEDITEC ELISA was calculated by adding 2 standard deviations to the mean OD of 20 replicate analyses of the negative control and was found to be 0.59 DU/mL (OD 450 nm = 0.034).

### 9.2 Diagnostic Specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte. It is 100%.

### 9.3 Diagnostic Sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte. It is 98%.

### 9.4 Method Comparison

The DEMEDITEC ELISA was compared with the Mikrogen Parvovirus B19 IgG 74 serum samples are assayed.

		Reference Assay	
		pos.	neg.
n = 74			
DEMEDITEC ELISA	pos.	55	0
	neg.	0	19

Agreement: 100%

### 9.5 Reproducibility

#### 9.5.1 Intra-assay

The intra-assay (within-run) precision of the DEMEDITEC ELISA was determined by 20x measurements of 12 samples covering the measuring range of the ELISA.

Sample	n	Mean OD	Intra-Assay CV (%)
1	20	0.03	9.7
2	20	0.26	5.7
3	20	0.07	6.0
4	20	0.99	8.8
5	20	0.76	5.9
6	20	0.71	9.9
7	20	1.23	6.8
8	20	1.37	3.7
9	20	1.25	3.4
10	20	1.54	3.9
11	20	2.00	3.6
12	20	2.03	1.7

#### 9.5.2 Inter-assay

The inter-assay variation of the DEMEDITEC ELISA was determined in 20 independent runs with 2 replicates per run.

Sample	n	Mean OD	Intra-Assay CV (%)
1	20	2.26	7.2
2	20	1.75	9.9
3	20	2.32	9.3

## 10 LIMITATIONS OF USE

Bacterial contamination or repeated freeze-thaw cycles of the specimen may affect the OD values. In immunocompromised patients and newborns serological data only have restricted value.

### 10.1 Interfering Substances

Haemoglobin (up to 4 mg/mL), Bilirubin (up to 0.5 mg/mL) and Triglyceride (up to 30 mg/mL) have no influence on the assay results.

## **11 LEGAL ASPECTS**

### **11.1 Reliability of Results**

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable national standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test. The test results are valid only if all controls are within the specified ranges and if all other test parameters are also within the given assay specifications. In case of any doubt or concern please contact DEMEDITEC.

### **11.2 Therapeutic Consequences**

Therapeutic consequences should never be based on laboratory results alone even if all test results are in agreement with the items as stated under point 11.1. Any laboratory result is only a part of the total clinical picture of a patient.

Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data. Only in cases where the laboratory results are in acceptable agreement with the overall clinical picture of the patient should therapeutic consequences be derived.

The test result itself should never be the sole determinant for deriving any therapeutic consequences.

### **11.3 Liability**

Any modification of the test kit and/or exchange or mixture of any components of different lots from one test kit to another could negatively affect the intended results and validity of the overall test. Such modification and/or exchanges invalidate any claim for replacement.

Claims submitted due to customer misinterpretation of laboratory results subject to point 11.2 are also invalid. Regardless, in the event of any claim, the manufacturer's liability is not to exceed the value of the test kit. Any damage caused to the test kit during transportation is not subject to the liability of the manufacturer.

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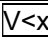
**SHORT INSTRUCTIONS FOR USE**

	All reagents and specimens must be allowed to come to room temperature (18 °C - 25 °C) before use.
	Dispense 100 µL of Controls into appropriate wells.
	Dispense 100 µL of sample into selected wells. <b>(Please note special sample treatment, point 5.3!)</b>
 60 min	Cover wells with foil. Incubate for <b>60 minutes</b> at 37 °C.
	Briskly shake out the contents of the wells.
	Rinse the wells <b>5 times</b> with diluted Wash Solution (300 µL per well).
	Strike the wells sharply on absorbent paper to remove residual droplets.
	Dispense 100 µL of Enzyme-Conjugate into each well.
 30 min	Incubate for <b>30 minutes</b> at room temperature.
	Briskly shake out the contents of the wells.
	Rinse the wells <b>5 times</b> with diluted Wash Solution (300 µL per well).
	Strike the wells sharply on absorbent paper to remove residual droplets.
	Add 100 µL of Substrate Solution to each well.
 15 min	Incubate for <b>15 minutes</b> at room temperature.
	Stop the reaction by adding 100 µL of Stop Solution to each well.
	Determine the optical density of each well at 450/620 nm.

## SYMBOLS USED WITH DEMEDITEC ASSAYS

Symbol	English	Deutsch	Française	Espanol	Italiano
	European Conformity	CE-Konformitätskennzeichnung	Conforme aux normes européennes	Conformidad europea	Conformità europea
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las Instrucciones	Consultare le istruzioni per l'uso
	In vitro diagnostic device	In-vitro-Diagnostikum	Diagnostic in vitro	Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Référence	Número de catálogo	No. di catalogo
	Lot. No. / Batch code	Chargen-Nr.	No. de lot	Número de lote	Lotto no
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Note warnings and precautions	Warnhinweise und Vorsichtsmaßnahmen beachten	Avertissements et mesures de précaution font attention	Tiene en cuenta advertencias y precauciones	Annoti avvisi e le precauzioni
	Storage Temperature	Lagerungstemperatur	Température de conservation	Temperatura de conservación	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeitsdatum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
<i>Distributed by</i>	Distributed by	Vertrieb durch	Distribution par	Distribución por	Distribuzione da parte di
	Version	Version	Version	Versión	Versione
	Single-use	Einmalverwendung	À usage unique	Uso único	Uso una volta